

Key factors that affect sonoporation efficiency in *in vitro* settings: The importance of standing wave in sonoporation

Manabu Kinoshita ^{a,*}, Kullervo Hynynen ^b

^a Department of Radiology, Brigham and Women's Hospital, and Harvard Medical School, 75 Francis Street, Boston, MA 02115, USA

^b Department of Medical Biophysics, University of Toronto and Imaging Research, Sunnybrook Health Sciences Centre, 2075 Bayview Avenue, Suite S6-65B, Toronto, Ont., Canada M4N 3M5

Received 17 May 2007

Available online 30 May 2007

Abstract

Ultrasound-induced intracellular drug delivery, sonoporation, is an appealing and promising technique for next generation drug delivery system. Many types of molecules, such as plasmid DNAs, siRNAs and peptides, have been demonstrated to be delivered into the cell by ultrasound with the aid of microbubbles both *in vitro* and *in vivo*. Although there are many reports on *in vitro* sonoporation, the efficiency of successful sonoporation and the viabilities of cells after the procedure documented in each report vary in a wide range, and the reasons for these differences are not fully understood. In this study, we have investigated how different experimental settings would affect sonoporation efficiency and cell viabilities after the procedure. Our results show that the fashion of cell culture (e.g. in suspension or in monolayer culture) and the presence of standing wave have a great impact on the overall results. These results indicate that *in vitro* sonoporation settings should be carefully evaluated in each experiment. The fact that standing wave is necessary to achieve high sonoporation efficiency may be a problematic issue for clinical application of sonoporation, as it may be difficult (although not impossible) to create standing wave in a human body.

© 2007 Elsevier Inc. All rights reserved.

Keywords: Sonoporation; Standing wave; siRNA; Gene therapy

Today the advancement of acoustic technology has made it possible to produce various kinds of bio-effects both *in vitro* and *in vivo* by ultrasound. These bio-effects include hyperthermia, such as coagulation of tissues [1–12], mechanical tissue destruction [13,14], disruption of tight junctions of the vasculature for enhanced local drug delivery [15–20], and sonoporation, a transient enhanced permeability of the cell membrane caused by the combination of ultrasound and microbubbles [21–38]. As all of these bio-effects can be created non-invasively *in vivo*, ultrasound is gathering a great attention as a next generation surgical or drug delivery system. Among them, sonopora-

tion has been vigorously investigated, as it allows plasmid DNAs, siRNAs, and peptides to be delivered intracellularly without the use of carriers or vectors both *in vitro* and *in vivo*. This concept has a great significance especially in the medical field, since it will widen the use of cell membrane impermeable agents. However, when we look into the published reports, sonoporation efficiency and the magnitude of cell toxicity, one of the most important factors for making sonoporation as a practical tool, varies in a very wide range. Some reports, including our previous reports, show only a small amount of successfully sonoporated cells accompanying a large number of non-viable cells [25,26,30,31,35]. Others report otherwise [34,38]. Although some reports have pursued the possibility of optimizing ultrasound parameters [26,35,39,40], there has been no report on investigating how the technical differences of experimental settings can impact the results. Solving this

* Corresponding author. Present address: Department of Neurosurgery, Osaka University Medical School, 2-2 Yamadaoka, Suita, Osaka 565-0871, Japan. Fax: +1 617 525 7450 (USA), +81 6 6879 3659 (Japan).

E-mail addresses: mkinosit@bwh.harvard.edu, m-kinoshita@nsurg.med.osaka-u.ac.jp (M. Kinoshita).

question is significantly important for refining sonoporation technology and for bringing it to practical application.

In this study, we have investigated how different experimental settings will affect sonoporation efficiencies and cell viabilities after the procedure. As reports on *in vitro* sonoporation experiment are mixed with cells sonicated in a suspension and cells sonicated in a monolayer culture, first we have investigated whether or not this difference will produce a difference in sonoporation efficiency. We have also tested if siRNA can be delivered by sonoporation in a monolayer cell culture. Finally, we investigated the role of standing wave in sonoporation. Standing wave is a condition where the reflected ultrasound beam from any sort of interface and the progressive ultrasound beam merge together, boosting the acoustic pressure in the field and altering the 3D structure of the ultrasound field from that created only by progressive ultrasound beam. Although producing standing wave in a human body is difficult, except for structures around the bone, in most of the reported sonoporation experiments, the presence or absence of standing wave is not clearly described. In addition, when standing wave is present, it is difficult to presume the true acoustic power applied to the sonicated samples. In this study, we have carefully created several ultrasound fields where standing wave is considered to be almost completely eliminated and evaluated its impact on sonoporation efficiency.

Materials and methods

Cell lines and culture. Rat C166 cells and C166-GFP cells, which are the wild-type and genetically manipulated cells stably expressing the enhanced green-fluorescent protein (EGFP), were obtained from the American Type Culture Collection (ATCC). C166 cells were maintained in Dulbecco's modified Eagle's medium (DMEM; Invitrogen, Carlsbad, CA) at 37 °C with 5% CO₂. Both media were supplemented with 10% fetal calf serum (FCS; Invitrogen). For C166-GFP cells, 0.2 mg/ml G418 (Geneticin, Invitrogen) was added in the medium.

Sonication device setup. Details of the experimental setup are shown in Fig. 1A. C166 or C166-GFP cells were cultured in a 6-well plate either in a suspension or in a monolayer fashion. Total number of cells in each well was fixed as 0.6×10^6 cells/1.5 ml serum-free DMEM. We used a 1.706-MHz planer ultrasound transducer with a surface area of 7 cm². The transducer's efficiency was equal to approximately 80% and it was used to convert the electrical power measured by the amplifier into acoustical power. The radiofrequency (RF) signals to drive the ultrasound transducers were generated by a pulsed-output ultrasound driving system (model UDS 04PF-CSA, Advanced Surgical Systems Inc., Tuscon, AZ). The transducer was submerged in degassed water and the samples were placed 6 cm above the surface of the transducer (Fig. 1A). The ultrasound parameters used were, 0–4 W/cm² acoustic power (AP) in continuous wave (CW; 100% duty cycle) mode for 15–120 s. Maximum temperature rise was less than 4.6 °C with 1.6 W/cm² AP \times CW 120 s. This was measured in 1.5 ml degassed water sonicated with standing present and microbubbles absent. The base line temperature was 36.5 °C. We also used OPTISON[®], a microbubble ultrasound contrast agent available from Amersham Health Inc. (Princeton, NJ), with a concentration of 5% or 10% (v/v) to aid cavitation production.

Sonoporation efficiency and cell viability assay for C166 cells. To determine the sonoporation efficiency of C166 cells under various ultrasound exposure conditions, cell membrane impermeable calcein (Sigma–Aldrich, Saint Louis, MO) was added to the cell suspension at

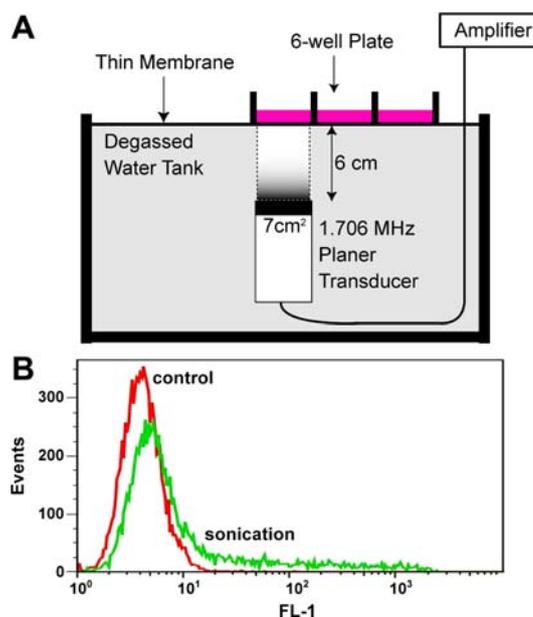


Fig. 1. Experimental setup for cell sonication using planer ultrasound. (A) A polystyrene 6-well plate was placed on a device designed for planer ultrasound exposure. Planer ultrasound (1.706 MHz center frequency) was delivered vertically to the cells from the bottom through a water tank and each sample was sonicated separately by moving the plate to the target of the focused ultrasound. (B) Sonoporated cells were detected using flow-cytometric analysis. Positive cells were labeled with calcein and detected on FL-1 channel.

a final concentration of 37.5 μ M. Cells were sonicated as mentioned above under various conditions and recovered, washed twice with phosphate-buffered saline (PBS) and analyzed by flow cytometry (FACScan; Becton–Dickinson, Mountain View, CA). Successfully sonoporated cells would uptake calcein and were detected by FL1 channel (Fig. 1B). To measure the viability of sonicated cells, we used CellTiter-Blue assay (Promega, Madison, WI). Sonicated cells were cultured for 2 h with the reagent and their viability was estimated by measuring the amount of reduced resorufin by its fluorescent intensity (ex: 560 nm, em: 590 nm) using a fluorescent microplate reader (SPECTRAMax GEMINI XS, Sunnyvale, CA).

siRNA transfection for EGFP suppression using microbubble-enhanced ultrasound. EGFP-specific siRNA; egfp siRNA (sense: 5'-GCA AGC UGA CCC UGA AGU UCA U-3', antisense: 5'-GAA CUU CAG GGU CAG CUU GCC G-3'), and control siRNA (sense: 5'-UUC UCC GAA CGU GUC ACG UdTdT-3', antisense: 5'-ACG UGA CAC GUU CGG AGA AdTdT-3') duplexes were purchased from Qiagen (Valencia, CA). The siRNA duplexes were suspended in the provided buffer solution and prepared according to the instruction supplied by the manufacturer. C166-GFP cells cultured in monolayer were suspended in their respective culture medium without serum. OPTISON[®] was added to the cell suspension at a final concentration of 5%. siRNAs were added at a concentration of 15 μ g/ml for gene silencing. After sonication, cells were recovered and 1.5 ml 10% serum and 50 U/ml penicillin (Invitrogen) supplemented medium was added and cultured for 48 h followed by flow-cytometric analysis of EGFP expression.

Flow-cytometric analysis. After cells were cultured for the indicated time, cells were harvested, recovered and washed twice with PBS. Flow-cytometric data collection was done using FACScan (Becton–Dickinson) and EGFP or calcein-positive cells were detected by FL1 channel. Acquired data was analyzed by FlowJo (Tree Star Inc, Ashland, OR). Two cell populations were identified on the forward scatter (FSC) vs. side scatter (SSC) plot and the population having larger values of both the FSC and SSC was gated as viable cells.

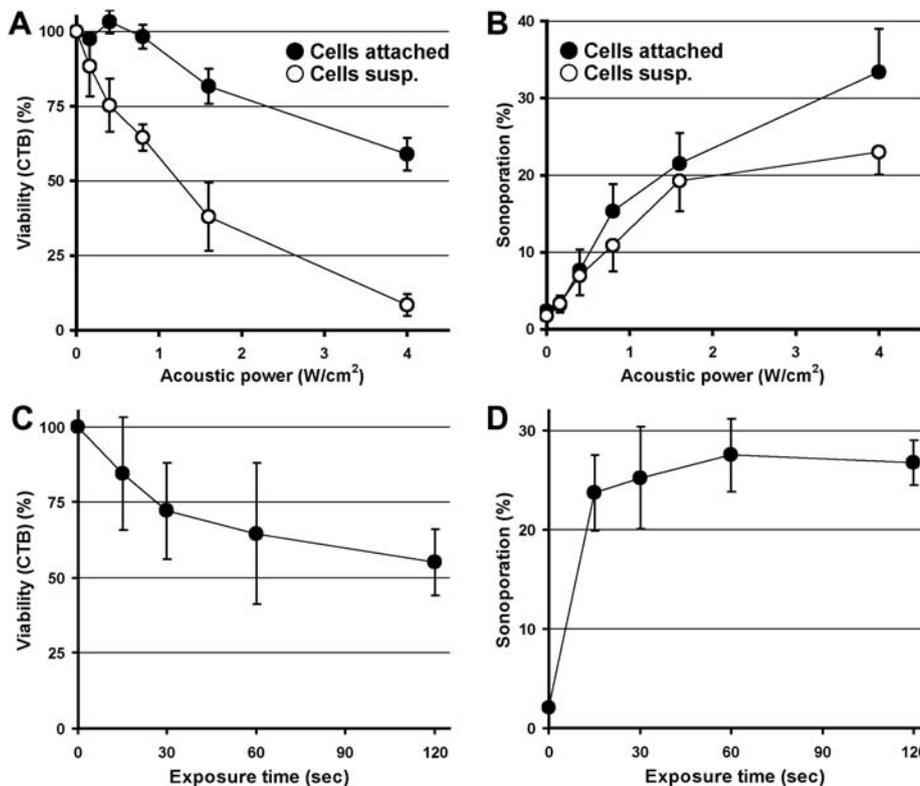


Fig. 2. Impact of the condition of cell culture on sonoporation. (A,B) Sonoporation with cells in suspension vs. attached to the dish. The concentration of OPTISON was 5% and the exposure time was fixed to 15 s. Cell viabilities (A) and sonoporation efficiencies among viable cells (B) are plotted as a function of the applied acoustic power. Data obtained from two conditions; cells in suspension (○) and cells cultured in monolayer (●) are presented. (C,D) Sonoporation efficiency as a function of exposure time. The concentration of OPTISON was increased up to 10% and the power was fixed to 1.6 W/cm² AP and the cells cultured in monolayer were sonoporated. Cell viabilities (C) and sonoporation efficiencies among viable cells (D) are plotted as a function of exposure time. All data are presented as means ± SD of three independent experiments.

Special settings to avoid standing wave production. For experiments where the production of standing wave had to be avoided, we built a water chamber above the cell culture samples with a thin plastic membrane separating the samples and the water chamber. Air bubbles in the cell culture medium were carefully removed, allowing the plastic membrane to have a direct contact with the medium. Three different settings were prepared as shown in Fig. 3.

Results and discussions

Although there are many reports on sonoporation [21–38], each experimental setting differs from one another and the reported sonoporation efficiency, which is one of the most significant information in sonoporation, varies. To understand the effect of the experimental settings themselves on sonoporation experiments, we first investigated the impact of the conditions of the cultured cells on sonoporation efficiency. When cells were sonicated in a suspension, although sonoporation was achieved among viable cells (Fig. 2B), majority of the cells were non-viable (Fig. 2A). On the other hand, when cells were cultured on the dish in a monolayer fashion, although the ratio of successfully sonoporated cells among viable cells did not differ from that in a condition with cells in suspension, the total viability was significantly higher than in suspensions (Fig. 2A and B). Increasing the concentration of

OPTISON from 5% to 10% did not contribute in improving the overall outcome (Fig. 2C), and extending the sonication time only decreased the total cell viability while sonoporation efficiency saturated at around 15–30 s CW exposure (Fig. 2C and D). This is a similar observation compared to what we have previously seen in cells sonicated in a suspension [30,31]. Delivery of siRNA against EGFP into cells with cell attached on the culture dish was also demonstrated possible (Fig. 4C) and approximately half of the cell sonoporated showed decreased EGFP production (Fig. 4C). These findings suggest that the conditions where cells are cultured have a great influence on the total outcome. As cells are structured in a matrix *in vivo*, it is assumed that sonoporation efficacy *in vivo* might be better compared to those performed *in vitro* with cell cultured in a suspension.

Next, we evaluated the impact of the presence or absence of standing wave on sonoporation efficiency. In this experiment, sonoporation was performed with cell attached on the bottom of the dish and standing wave was eliminated by building a water chamber above the sonication samples, allowing the progressive ultrasound wave to propagate through the sonication samples without causing reflected ultrasound wave (Fig. 3, left). When standing wave was nearly completely eliminated from the

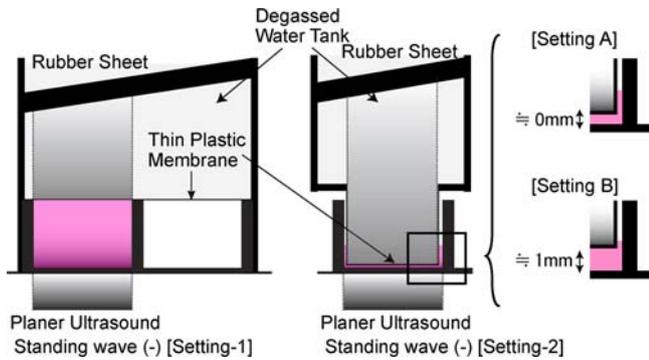


Fig. 3. Experimental setting eliminating standing wave. Experimental settings eliminating standing wave. Cells were cultured in a monolayer fashion on the bottom of the dish. A water chamber was constructed above the samples which allow the ultrasound beam to propagate through the sample without causing standing wave. In some cases, the interface of the water chamber and the sonicated cells was settled close to 0 mm (setting A) or 1 mm above (setting B) the bottom of the dish.

experimental settings, only a negligible amount of cells achieved sonoporation, even under condition where maximum sonoporation was achieved in the presence of standing wave (Fig. 4B). Cell viability after sonoporation was nearly 100% even after 5 W/cm^2 exposures (Fig. 4A). As there is a possibility that all the OPTISON[®] microbubbles were pushed away from the cultured cells during sonication, we prepared two different settings in addition. One with the water chamber directly contacting the cells attached on the bottom of the wells (Fig. 3, right setting-2A), and the other lifting the water chamber 1 mm above the cells attached on the bottom (Fig. 3, right setting-2B). Again, when standing wave was absent, both of the settings failed to achieve a satisfactory amount of sonoporation, even up to 5 W/cm^2 exposures (Fig. 4A and B). These results suggest that creation of standing wave is one of the key elements in achieving successful sonoporation in

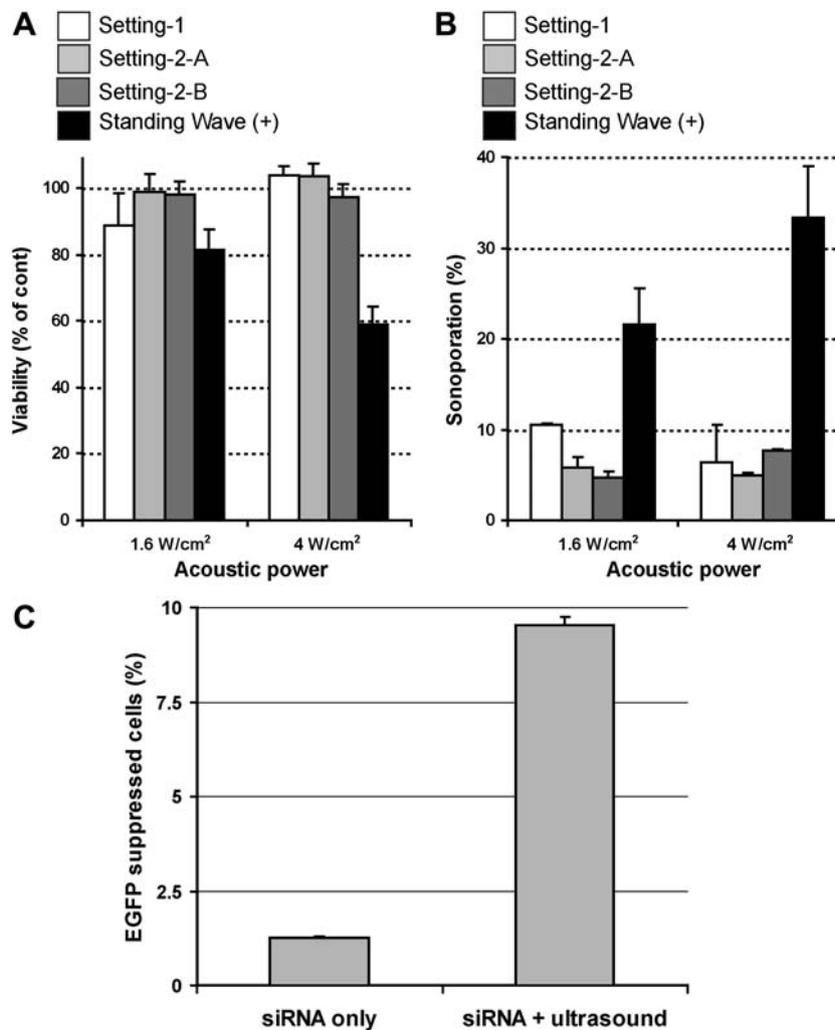


Fig. 4. The impact of standing wave on sonoporation and inhibition of EGFP stable expression in C166-GFP cells by microbubble-enhanced ultrasound delivered siRNA. (A,B) Sonoporation without the presence of standing wave. The concentration of OPTISON was 5% and the exposure time was fixed to 15 s. Note that both loss of cell viability (A) and sonoporation efficacy (B) have dramatically dropped when standing wave is absent. All data are presented as means \pm SD of three independent experiments. (C) EGFP gene suppression with intracellular siRNA delivery. The conditions were 5% OPTISON with 1.6 W/cm^2 AP \times 15 s exposure and the cells were cultured in a monolayer fashion. The concentration of siRNA was the same as previously reported where cells were sonicated in suspension [30]. Approximately half of the sonoporated cells underwent gene suppression (means \pm SD of three independent experiments).

in vitro experiments. This finding could be crucial in applying sonoporation to real clinical use, as the production of standing wave may be difficult in the human body and special arrangements may be needed to establish the desired acoustic fields. When we look into the reported sonoporation experiments *in vivo*, the animals often used for demonstration are rats or rabbits [27,29,33,34,37]. Again, the animals are sonicated under a condition where it is difficult to assess the true ultrasound field and the presence or absence of standing wave is not clear. However, based on the data available, it is likely that standing waves were generated in many of these experiments. As the information on the relationship between the applied ultrasound field characteristics and the achieved sonoporation efficiency is crucial to render this technique practical in the future clinical application, we suggest that researchers should be aware of the above mentioned matters in performing and reporting sonoporation experiments.

Acknowledgments

This investigation was supported by the Shinya International Exchange Fund, the Osaka Medical Research Foundation for Incurable Diseases, the Osaka Neurological Institute and NIH Grant (EB003268).

References

- [1] J.Y. Chapelon, J. Margonari, F. Vernier, F. Gorry, R. Ecochard, A. Gelet, In vivo effects of high-intensity ultrasound on prostatic adenocarcinoma Dunning R3327, *Cancer Res.* 52 (1992) 6353–6357.
- [2] R. Chopra, M. Burtnyk, M.A. Haider, M.J. Bronskill, Method for MRI-guided conformal thermal therapy of prostate with planar transurethral ultrasound heating applicators, *Phys. Med. Biol.* 50 (2005) 4957–4975.
- [3] C.J. Diederich, Thermal ablation and high-temperature thermal therapy: overview of technology and clinical implementation, *Int. J. Hyperthermia* 21 (2005) 745–753.
- [4] W.J. Fry, J.W. Barnard, F.J. Fry, J.F. Brennan, Ultrasonically produced localized selective lesions in the central nervous system, *Am. J. Phys. Med.* 34 (1955) 413–423.
- [5] J.D. Hazle, C.J. Diederich, M. Kangasniemi, R.E. Price, L.E. Olsson, R.J. Stafford, MRI-guided thermal therapy of transplanted tumors in the canine prostate using a directional transurethral ultrasound applicator, *J. Magn. Reson. Imaging* 15 (2002) 409–417.
- [6] N.T. Sanghvi, R.S. Foster, R. Bihrl, R. Casey, T. Uchida, M.H. Phillips, J. Syrus, A.V. Zaitsev, K.W. Marich, F.J. Fry, Noninvasive surgery of prostate tissue by high intensity focused ultrasound: an updated report, *Eur. J. Ultrasound* 9 (1999) 19–29.
- [7] G. ter Haar, I. Rivens, L. Chen, S. Riddler, High intensity focused ultrasound for the treatment of rat tumours, *Phys. Med. Biol.* 36 (1991) 1495–1501.
- [8] F. Wu, Z.B. Wang, H. Zhu, W.Z. Chen, J.Z. Zou, J. Bai, K.Q. Li, C.B. Jin, F.L. Xie, H.B. Su, Feasibility of US-guided high-intensity focused ultrasound treatment in patients with advanced pancreatic cancer: initial experience, *Radiology* 236 (2005) 1034–1040.
- [9] N. McDannold, M. Moss, R. Killiany, D.L. Rosene, R.L. King, F.A. Jolesz, K. Hynynen, MRI-guided focused ultrasound surgery in the brain: tests in a primate model, *Magn. Reson. Med.* 49 (2003) 1188–1191.
- [10] C.M. Tempany, E.A. Stewart, N. McDannold, B.J. Quade, F.A. Jolesz, K. Hynynen, MR imaging-guided focused ultrasound surgery of uterine leiomyomas: a feasibility study, *Radiology* 226 (2003) 897–905.
- [11] K. Hynynen, O. Pomeroy, D.N. Smith, P.E. Huber, N.J. McDannold, J. Kettenbach, J. Baum, S. Singer, F.A. Jolesz, MR imaging-guided focused ultrasound surgery of fibroadenomas in the breast: a feasibility study, *Radiology* 219 (2001) 176–185.
- [12] K. Hynynen, G.T. Clement, N. McDannold, N. Vykhodtseva, R. King, P.J. White, S. Vitek, F.A. Jolesz, 500-element ultrasound phased array system for noninvasive focal surgery of the brain: a preliminary rabbit study with ex vivo human skulls, *Magn. Reson. Med.* 52 (2004) 100–107.
- [13] N.B. Smith, K. Hynynen, The feasibility of using focused ultrasound for transmyocardial revascularization, *Ultrasound Med. Biol.* 24 (1998) 1045–1054.
- [14] B.C. Tran, J. Seo, T.L. Hall, J.B. Fowlkes, C.A. Cain, Microbubble-enhanced cavitation for noninvasive ultrasound surgery, *IEEE Trans. Ultrason. Ferroelectr. Freq. Control* 50 (2003) 1296–1304.
- [15] K. Hynynen, N. McDannold, N. Vykhodtseva, F.A. Jolesz, Noninvasive MR imaging-guided focal opening of the blood–brain barrier in rabbits, *Radiology* 220 (2001) 640–646.
- [16] K. Hynynen, N. McDannold, N.A. Sheikov, F.A. Jolesz, N. Vykhodtseva, Local and reversible blood–brain barrier disruption by noninvasive focused ultrasound at frequencies suitable for trans-skull sonications, *Neuroimage* 24 (2005) 12–20.
- [17] K. Hynynen, N. McDannold, N. Vykhodtseva, S. Raymond, R. Weissleder, F.A. Jolesz, N. Sheikov, Focal disruption of the blood–brain barrier due to 260-kHz ultrasound bursts: a method for molecular imaging and targeted drug delivery, *J. Neurosurg.* 105 (2006) 445–454.
- [18] N. McDannold, N. Vykhodtseva, S. Raymond, F.A. Jolesz, K. Hynynen, MRI-guided targeted blood-brain barrier disruption with focused ultrasound: histological findings in rabbits, *Ultrasound Med. Biol.* 31 (2005) 1527–1537.
- [19] M. Kinoshita, N. McDannold, F.A. Jolesz, K. Hynynen, Targeted delivery of antibodies through the blood-brain barrier by MRI-guided focused ultrasound, *Biochem. Biophys. Res. Commun.* 340 (2006) 1085–1090.
- [20] M. Kinoshita, N. McDannold, F.A. Jolesz, K. Hynynen, Noninvasive localized delivery of Herceptin to the mouse brain by MRI-guided focused ultrasound-induced blood–brain barrier disruption, *Proc. Natl. Acad. Sci. USA* 103 (2006) 11719–11723.
- [21] K. Anwer, G. Kao, B. Proctor, I. Anscombe, V. Florack, R. Earls, E. Wilson, T. McCreery, E. Unger, A. Rolland, S.M. Sullivan, Ultrasound enhancement of cationic lipid-mediated gene transfer to primary tumors following systemic administration, *Gene Ther.* 7 (2000) 1833–1839.
- [22] D.L. Miller, S.V. Pislaru, J.E. Greenleaf, Sonoporation: mechanical DNA delivery by ultrasonic cavitation, *Somat. Cell Mol. Genet.* 27 (2002) 115–134.
- [23] H. Azuma, N. Tomita, Y. Kaneda, H. Koike, T. Ogihara, Y. Katsuoka, R. Morishita, Transfection of NF κ B-decoy oligodeoxynucleotides using efficient ultrasound-mediated gene transfer into donor kidneys prolonged survival of rat renal allografts, *Gene Ther.* 10 (2003) 415–425.
- [24] C.X. Deng, F. Sieling, H. Pan, J. Cui, Ultrasound-induced cell membrane porosity, *Ultrasound Med. Biol.* 30 (2004) 519–526.
- [25] H.R. Guzman, D.X. Nguyen, A.J. McNamara, M.R. Prausnitz, Equilibrium loading of cells with macromolecules by ultrasound: effects of molecular size and acoustic energy, *J. Pharm. Sci.* 91 (2002) 1693–1701.
- [26] H.R. Guzman, A.J. McNamara, D.X. Nguyen, M.R. Prausnitz, Bioeffects caused by changes in acoustic cavitation bubble density and cell concentration: a unified explanation based on cell-to-bubble ratio and blast radius, *Ultrasound Med. Biol.* 29 (2003) 1211–1222.
- [27] N. Hashiya, M. Aoki, K. Tachibana, Y. Taniyama, K. Yamasaki, K. Hiraoka, H. Makino, K. Yasufumi, T. Ogihara, R. Morishita, Local delivery of E2F decoy oligodeoxynucleotides using ultrasound with microbubble agent (Optison) inhibits intimal hyperplasia after

- balloon injury in rat carotid artery model, *Biochem. Biophys. Res. Commun.* 317 (2004) 508–514.
- [28] H. Inagaki, J. Suzuki, M. Ogawa, Y. Taniyama, R. Morishita, M. Isobe, Ultrasound-microbubble-mediated NF-kappaB decoy transfection attenuates neointimal formation after arterial injury in mice, *J. Vasc. Res.* 43 (2006) 12–18.
- [29] P.E. Huber, M.J. Mann, L.G. Melo, A. Ehsan, D. Kong, L. Zhang, M. Rezvani, P. Peschke, F. Jolesz, V.J. Dzau, K. Hynynen, Focused ultrasound (HIFU) induces localized enhancement of reporter gene expression in rabbit carotid artery, *Gene Ther.* 10 (2003) 1600–1607.
- [30] M. Kinoshita, K. Hynynen, A novel method for the intracellular delivery of siRNA using microbubble-enhanced focused ultrasound, *Biochem. Biophys. Res. Commun.* 335 (2005) 393–399.
- [31] M. Kinoshita, K. Hynynen, Intracellular delivery of Bak BH3 peptide by microbubble-enhanced ultrasound, *Pharm. Res.* 22 (2005) 716–720.
- [32] M. Shimamura, N. Sato, Y. Taniyama, S. Yamamoto, M. Endoh, H. Kurinami, M. Aoki, T. Ogihara, Y. Kaneda, R. Morishita, Development of efficient plasmid DNA transfer into adult rat central nervous system using microbubble-enhanced ultrasound, *Gene Ther.* 11 (2004) 1532–1539.
- [33] M. Shimamura, N. Sato, Y. Taniyama, H. Kurinami, H. Tanaka, T. Takami, T. Ogihara, M. Tohyama, Y. Kaneda, R. Morishita, Gene transfer into adult rat spinal cord using naked plasmid DNA and ultrasound microbubbles, *J. Gene. Med.* 7 (2005) 1468–1474.
- [34] S. Sonoda, K. Tachibana, E. Uchino, A. Okubo, M. Yamamoto, K. Sakoda, T. Hisatomi, K.H. Sonoda, Y. Negishi, Y. Izumi, S. Takao, T. Sakamoto, Gene transfer to corneal epithelium and keratocytes mediated by ultrasound with microbubbles, *Invest. Ophthalmol. Vis. Sci.* 47 (2006) 558–564.
- [35] V.G. Zarnitsyn, M.R. Prausnitz, Physical parameters influencing optimization of ultrasound-mediated DNA transfection, *Ultrasound Med. Biol.* 30 (2004) 527–538.
- [36] S. Tsunoda, O. Mazda, Y. Oda, Y. Iida, S. Akabame, T. Kishida, M. Shin-Ya, H. Asada, S. Gojo, J. Imanishi, H. Matsubara, T. Yoshikawa, Sonoporation using microbubble BR14 promotes pDNA/siRNA transduction to murine heart, *Biochem. Biophys. Res. Commun.* 336 (2005) 118–127.
- [37] Y. Taniyama, K. Tachibana, K. Hiraoka, T. Namba, K. Yamasaki, N. Hashiya, M. Aoki, T. Ogihara, K. Yasufumi, R. Morishita, Local delivery of plasmid DNA into rat carotid artery using ultrasound, *Circulation* 105 (2002) 1233–1239.
- [38] M. Duvshani-Eshet, L. Baruch, E. Kesselman, E. Shimoni, M. Machluf, Therapeutic ultrasound-mediated DNA to cell and nucleus: bioeffects revealed by confocal and atomic force microscopy, *Gene Ther.* 13 (2006) 163–172.
- [39] A.A. Rahim, S.L. Taylor, N.L. Bush, G.R. ter Haar, J.C. Bamber, C.D. Porter, Spatial and acoustic pressure dependence of microbubble-mediated gene delivery targeted using focused ultrasound, *J. Gene. Med.* 8 (2006) 1347–1357.
- [40] A. Rahim, S.L. Taylor, N.L. Bush, G.R. ter Haar, J.C. Bamber, C.D. Porter, Physical parameters affecting ultrasound/microbubble-mediated gene delivery efficiency in vitro, *Ultrasound Med. Biol.* 32 (2006) 1269–1279.